V E N O M

BY **HYPERRICE**

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
**WARNINGS AND CAUTION**

- Do not immerse in water. Keep liquids away from unit. Clean only with a damp cloth. Do not disassemble.
- Turn off immediately and remove from body at first indication of discomfort.
- Do not operate continuously for more than 1 hour. After 1 hour use, allow unit to rest for 30 min. before re-using.
- Unplug the charger from battery after charging or prior to use. Venom will not function while charger is plugged into battery.
- For storage or travel, remove battery and store in supplied carton box.
- The Venom is to be used only as intended by adults or children over 8 years old with adult supervision and only after reading the supplied instructions for use and warnings.
- The Venom is designed to get hot during normal operation. Persons sensitive to heat must exercise caution when using.
- The Venom must only be used with the voltage indicated by the supplied battery (Model 42000-001-00) and the approved supplied charger indicated with "V" logo.

DISCLAIMER

PLEASE DO NOT USE THE VENOM, OR ANY HEAT AND VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals when using this device. These contra indications do not mean that you are not able to use the Venom but we do advise you to consult a doctor first. Ongoing research is being done to the effects of heat and vibration for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of heat and vibration into a treatment plan is advisory in a number of cases. This must be done on the advice of a doctor, specialist or physiotherapist.

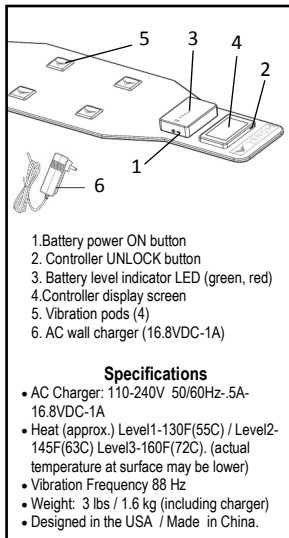1. Battery power ON button
2. Controller UNLOCK button
3. Battery level indicator LED (green, red)
4. Controller display screen
5. Vibration pods (4)
6. AC wall charger (16.8VDC-1A)

Specifications

- AC Charger: 110-240V 50/60Hz-.5A-16.8VDC-1A
- Heat (approx.) Level1-130F(55C) / Level2-145F(63C) Level3-160F(72C). (actual temperature at surface may be lower)
- Vibration Frequency 88 Hz
- Weight: 3 lbs / 1.6 kg (including charger)
- Designed in the USA / Made in China.

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